



BON CHARGE

**RED LIGHT  
TOOTHBRUSH™**

## Indications For Use

---

- For daily maintenance of oral hygiene
- To support healthy teeth and gums as part of a regular brushing routine
- To gently stimulate and massage gums with near-infrared and red light and sonic vibration.

## Device Description

---

The BON CHARGE Red Light Toothbrush combines high-frequency sonic vibration with red light (660 nm) and near-infrared (830 nm) to enhance oral hygiene and support gum care.

Depending on the selected mode, the toothbrush is engineered to deliver sonic vibration with the benefits of near-infrared and red light, offering effective cleaning and gentle gum stimulation.

## Important Safety Information





---

- Keep out of reach of children.
- Do not chew or bite the brush.
- Use only for cleaning teeth, gums, and tongue.
- Avoid shining light directly into the eyes.
- Consult a healthcare professional if you experience discomfort.
- Fully charge the toothbrush before first use.
- Do not use the device or accessories if damaged.
- Consult your dental professional before use if you have oral health concerns.
- Each family member should use their own toothbrush and brush head.

# Contents

QTY	PRODUCT
1 x	Toothbrush Handle
2 x	Replacement Brush Heads
1 x	Charging Base
1 x	USB Charging Cable
1 x	QR Code (Link to Online User Manual)
1 x	User Manual

## Product Features


LIGHT MODES	WAVELENGTHS	LED TECHNOLOGY
3 - Red light 3 - Near-infrared	Red Light: 660nm Near-Infrared: 830nm	2 Wavelengths 6 LEDs total
BRUSHING MODES	CLEAN (320HZ) 	SENSITIVE (310HZ) 
4 Modes	Powerful sonic vibration with near-infrared and red light.	Sensitive sonic clean with near-infrared and red light.
	CARE (315Hz) 	POLISH (290HZ) 
	Gentle sonic clean with near-infrared and red light.	Sonic clean with polishing action, near-infrared, and red light.

## MEMORY FUNCTION

The toothbrush recalls the last-used mode on startup

## INTELLIGENT BRUSHING

- 2-minute smart timer (Dentist recommended)
- 30-second interval pause to prompt you to switch to different areas of your mouth

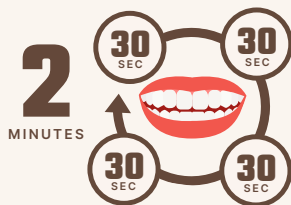

## Battery & Charging

### LONG-LASTING

- Fully charges in ~10 hours
- Lasts up to 4 weeks with twice-daily 2-minute use


### CHARGING

1. Ensure the toothbrush is turned off.
2. Plug the USB charging cable into the charging base and a USB power source\* (not included).
3. Place the toothbrush handle onto the charging base.

### CHARGING INDICATOR

The white charging indicator will flash while charging and remain solid when the battery is fully charged.

- Charging: **flashing WHITE.**
- Charged: **solid WHITE.**

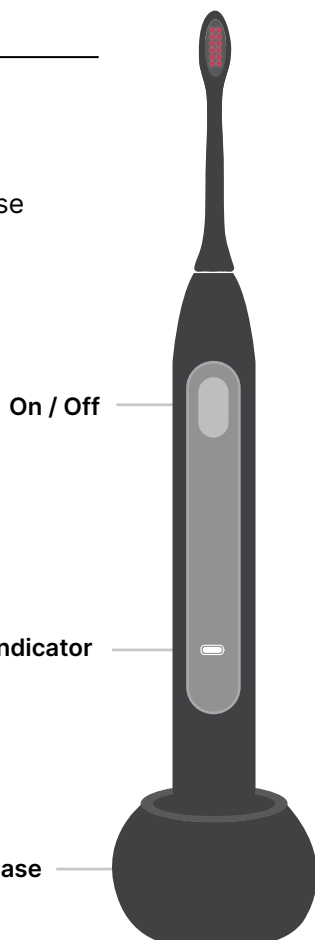

## CHARGING TIPS

- Charge fully before first use
- Fully recharge if not used for 3+ months or the battery is drained
- To maintain battery life, store on the charging base when not in use
- It's normal for the handle to feel warm while charging

***\*IMPORTANT:** This device is designed for use with a standard USB charger (DC 5V/1A). Fast-charging USB adapters may damage the product. Only use the USB cable supplied with the product.*

## How To Use

STEP 1:	STEP 2:
<b>Attach the brush head</b> to the handle. A small gap will remain between the head and the body of the toothbrush.	<b>Apply toothpaste</b> and wet the bristles.
STEP 3:	STEP 4:
<b>Insert the brush</b> into your mouth before switching on to avoid splashing.	<b>Press the power button</b> to cycle through modes.
STEP 5:	STEP 6:
Hold the brush at a <b>45° angle to the gumline</b> (especially in Care and Polish modes). Allow for gentle contact between the brush and your teeth.	<b>Gently move from tooth to tooth.</b> The timer will pause every 30 seconds to prompt you to move to the next quadrant of your mouth.

***Note:** Some users may experience slight gum bleeding during initial use. This usually resolves in a few days. If bleeding continues after two weeks, consult your dentist or healthcare professional.*

## Cleaning & Storage

---

- Rinse the brush head thoroughly after use
- Remove the brush head and rinse the shaft of the handle
- Do not cover the head after cleaning—allow to air dry
- Store upright in a cool, dry place, away from sunlight
- Unplug if not in use for an extended period

### REPLACEABLE LED ELECTRIC TOOTHBRUSH HEADS

The replaceable toothbrush heads are designed exclusively for use with the BON CHARGE Red Light Toothbrush. Replace every 3 months or sooner if worn, discolored, or if you notice a change in performance.

#### TO REMOVE THE HEAD:

- Hold the handle firmly.
- Pull the brush head straight off. Do not twist or bend.

#### CARE TIPS:

- Rinse with warm (not hot) water
- Avoid applying excessive pressure while brushing
- Let brush heads dry fully between uses

To reorder, visit: [boncharge.com](https://boncharge.com)

*The American Dental Association recommends replacing toothbrush heads after three months of use or sooner, especially if you have been sick or if you have stored it near toothbrushes of family members who have been sick.*

# Parameters

---

Rated Voltage:	DC 3.7V
Input Voltage:	DC 5V/1A
Standby Power:	≤0.5W
Battery:	800mAh Li-ion
Battery Life:	4+ weeks per charge
Vibration Frequency:	290 Hz - 320 Hz
Waterproof Level:	IPX7

# Warnings

---

- **STOP** use and consult a healthcare professional if you experience discomfort, irritation, dizziness, or other adverse reactions
- **DO NOT** stare at light source or direct light onto the eyes. This may cause temporary loss of sight.
- **DO NOT** use over open wounds or infected areas
- **DO NOT** use if device is damaged, dropped, or malfunctioning
- **ONLY** use accessories supplied or approved by the manufacturer
- **DO NOT** attempt to repair, modify, or disassemble—no serviceable or replaceable parts inside
- **DO NOT** submerge in water or charge near water
- **DO NOT** use if you've had recent oral surgery or are recovering from infection, as this may result in the worsening of symptoms.  
**DO NOT** use the device in any other way than specified in this manual.

## Danger (Electric Shock Risk)

---

- If the device falls into water while charging, unplug it before retrieving
- **DO NOT** use or charge near water
- **DO NOT** handle charger with wet hands
- Always unplug by pulling the plug, not the cable
- **DO NOT** modify, tamper with, or disassemble the charger or device
- **DO NOT** modify the charger

## Contraindications

---

**DO NOT** use the device if:

- You have active cancer or a cancerous lesion
- You have known allergies to light, plastic, or metal
- You are taking medication that causes photosensitivity
- You have a light-induced headache history
- You suffer from epilepsy or seizures triggered by light
- You suffer from any photosensitive disorder (sensitization to light)

## Light Sensitivity Test

---

Before using your toothbrush for the first time, we recommend conducting a sensitivity test to check for any potential issues.

1. Place device just above your inner forearm
2. Turn on and leave in position for **10 minutes**.
3. If discomfort occurs, stop use immediately
4. Observe the skin over the next 24 hours
  - a. If redness or irritation occurs, do not use the device
  - b. Consult your healthcare professional



# Troubleshooting

---

PROBLEM	POSSIBLE CAUSE	SOLUTION
The toothbrush will not work	May require charging	Check that the toothbrush is charged. If not, fully charge the toothbrush.
The toothbrush is charged but will not start	Toothbrush head not in place	Check that the toothbrush head is inserted all the way onto the toothbrush leaving only a small gap. If not, gently remove and re-insert the toothbrush head onto the toothbrush and try to switch it on again.
The toothbrush is not charging	The toothbrush may not be properly inserted into the charging base. The charging base may not be connected to a power source or cables may be damaged/ require replacement.	Check that the toothbrush is fully inserted onto the charging base and that the charging base is connected to a power source. To check, safely disconnect from the power source and inspect cables for damage. If damaged, do not operate the unit and contact Customer Support.
The light is not staying on while brushing	Toothbrush head not in place/poor contact	Check that the toothbrush head is inserted all the way onto the toothbrush. Cycle between modes to check that there is power/ contact between the toothbrush and toothbrush head in place. If lights are still not coming on after 10 seconds of cycling between

<b>The brush is not cycling between modes</b>	Toothbrush button may not have been pressed and continuously held down to reach the preferred mode	Check that you are pressing the main toothbrush button in sequence to toggle between the 4 different cleaning modes. Switch it off. Next time you switch on, the toothbrush will start from the last used cleaning mode.
<b>Increased pressure is required to achieve sufficient cleaning</b>	Toothbrush head may be worn or damaged	Check for any signs of toothpaste buildup, damage, bending of the bristles or scratch marks to the toothbrush head. Replace the brush head at least every three months or when worn or damaged.

## Disposal

---

Do not dispose of this product in household waste. It must be taken to a designated e-waste collection point. Inform staff the device contains a lithium battery.

## Warranty

---

- 1-year limited warranty from date of purchase
- Covers defects due to faulty materials or workmanship
- Proof of purchase required

Warranty **void** if:

Brush heads are damaged due to misuse or alteration

Product shows wear and tear (scratches, discoloration, etc.)

Damage caused disassembling the product, negligence, or caused by using an attachment not intended for use with this product.

To obtain service within the warranty period, contact Customer Support at **hello@boncharge.com**. Do not return to the address on the original packaging.

## Disclaimer

---

**This product is not intended to diagnose, cure or prevent any disease.**


















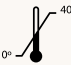


BON CHARGE makes no claims, representations or warranties regarding the ability of this product to cure any physical ailments. A qualified Healthcare Professional should be consulted with regard to any condition requiring medical attention.

## Customer Support

---

For questions or tips on troubleshooting, brushing modes, care tips, warranties or returns, please contact our support team at **hello@boncharge.com**.

# Labels & Symbols

SYMBOLS	DESCRIPTION	SYMBOLS	DESCRIPTION
	Store away from sunlight		Manufacturer
	Meets safety and performance standards for sale in Europe		Complies with FCC limits for electromagnetic emissions
	Do not stare at light source		Consult instructions for use
	Manufacture date		Recycle cardboard
	Safe to use and store only within these temperature ranges		Use and store the device only within these humidity
	Keep Dry		Product distributor
	Do not dispose in household waste. Separate collection for electrical and electronic equipment.		RoHS directive 2011/65/ EU - Restriction of the use of certain hazardous substances in electrical and electronic equipment
	UK Conformity Assessed as per UK Directives on health, safety and environmental protection standards		Lithium-ion battery is part of the device
	Only use USB cable provided		Unique Device Identifier
	Lot number		EU Representative



**Light Tree Ventures Europe B.V.**  
Laan Van Ypenburg 108, 2497 GC,  
The Hague, The Netherlands



**BLUBLOX PTY LTD**  
2/40 King Edward Road,  
Osborne Park WA 6017

**Contact:** [hello@boncharge.com](mailto:hello@boncharge.com)

**Website:** [www.boncharge.com](http://www.boncharge.com)